

## BEDPAN LINER GRADUATED 2000ML

**Eco-Friendly Medical Pulp** 





## Bedpan Liner Graduated 2000ml Eco-Friendly Medical Pulp







The Bedpan Liner Graduated with 2L rated capacity is made from eco-friendly biodegradable kraft material, intended for use as receptacle for urine and/or faeces. To use the bedpan liner, place it over a plastic or moulded pulp support. Alternatively, the bedpan liner can also be placed under commode chair for toileting.

Manufactured as single-use disposable medical pulp to meet the requirement to hold water and/or human waste for a minimum period of four hours at an initial temperature of 35°C. Suitable for disposal via macerator which complies with PAS29:1999 "Disposable Pulp For Use in Healthcare".

## Suitable for







**Nursing Homes** 



Personal Use




**Medical Clinics** 



**Disability Homes** 


## **User Guide**

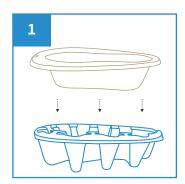


Place Bedpan Liner over the plastic or disposable moulded pulp bedpan support.



Slide the Bedpan Liner with support under the buttock. To prevent spills, place a sheet of incontinence pad under the bedpan.

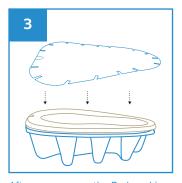

After use, remove the Bedpan Liner with support, place a lid cover over it to contain the odour and sight.

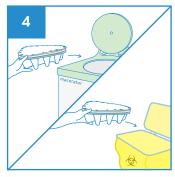

Dispose the Bedpan Liner with disposable support directly into macerator or empty the content from Bedpan Liner before disposing into appropriate waste bins.

Note: DO NOT dispose the plastic support into macerator.



Stackable



Water Retention



Macerator Friendly



**Eco-Friendly** 



Recycled

PRODUCT DESCRIPTION

**PRODUCT SIZE** 

**PACKAGING** 

**CARTON DIMENSION** 

**ORDER CODE** 

Bedpan Liner Graduated 2000ml

 $360 \times 300 \times 100 \text{ mm}$ 

100 pcs per carton

755 x 370 x 310 mm

100APPBPG2000







Manufactured in accordance with PAS29:1999
Disposable Pulp Products For Use in Healthcare

